

DAKEWE



Sureprint S200 Slide Printer

Dual input, dual output



SurePrint S200 slide printer




Featuring non-contact UV laser marking technology, Dakewe SurePrint S200 slide printer offers fast, rich-content and consumable-free printing. The innovative upward slide delivery system allows smooth and safe transport of slides. The printing is crisp and permanent, resists processing solvents, heat and scratch, reducing errors and enabling traceability.

High tech printing technology.

Thanks to the non-contact UV laser marking technology, S200 offers:

- High printing speed – as fast as 5 seconds per slide
- Crisp print – 2500 dpi high resolution and great scan rate
- Minimal maintenance – non-contact printing results in little to none maintenance
- Rich content – Alphanumeric, multiple languages, symbols, 1D & 2D codes, graphics

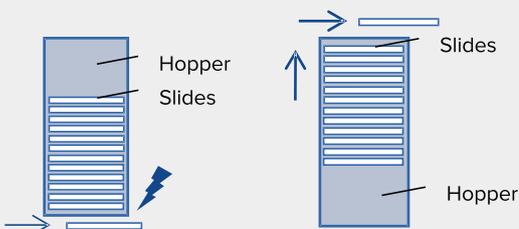

Left: downward slide delivery
Right: upward slide delivery in S200

Innovative upward slide delivery.

Gravity-driven delivery system are prone to break the slides and induce blockage.

This will never happen in S200 as it implants upward slide delivery system, slides are delivered upwards and then horizontally to the printing area, reducing the possibility of slide fragmentation.

Dual input, dual output.

- **Dual input, smart selection** – with two loading hoppers, S200 can store two types of slides separately and automatically print the requested slide in the corresponding hopper
- **Large capacity** – 200-slide loading capacity offers relief from repetitive reloading
- **Dual output, double printing modes** – upper slot for on-demand printing and lower hopper for batch printing, accommodating your laboratory routine



Fit seamlessly into your current workflow

- **Slide compatibility** – compatible with most of slides with paint coated tab
- **Seamless integration** – send printing command from the built-in software, PC software or integrated LIS
- **Easy operation** – windows-based system with concise UI

Easy and safe operation.

- **User safety** – triple-layer filters in the air purification system keep odor and dust away from operator
- **Front loading feature** – convenient loading and unloading



Specifications

Printing technology	Non-contact UV laser marking technology
Print mode	On-demand and batch printing
Loading hoppers	2, up to 100 slides in each hopper
Output hopper	Up to 50 slides
Print speed	Up to 5 seconds/slide
Print options	Alphanumeric, multiple languages, symbols, 1D & 2D codes, graphics
Air purification	Built-in activated carbon air purification system
System interface	LIS, HIS, PIS, PACS
Printer interface	USB and Ethernet
Power requirements	AC 100 - 240 V, 50/60 Hz, 210 VA
Dimensions (W×D×H)	240 × 440 × 450 mm
Weight	Approx. 32.5 kg
Operating systems	Win 7, Win 10, Win 11, Linux, macOS



DAKEWE

Dakewe Medical
Shenzhen, China
www.dakewemedical.com
enquiry@dakewe.com

SUREPRINT S200 SLIDE PRINTER
Ver 1.3, Rev. Apr 2024

